

# CR-X

**RMC**  
**Products**  
by Boeckeler®

## The Cost-Effective Solution to Cryosectioning

Vibrationless LN<sub>2</sub>  
Delivery

Side fill for easy  
access to chamber

Four Temperature Presets



Cryo Workstation  
(not shown) for easy  
section manipulation


Four knife holders  
included for flexible  
combination of glass and  
diamond knives  
and trimming tools

Mounts and dismounts on almost any  
ultramicrotome in seconds:  
RMC Models: PT-XL, PT-X, PT-PC  
(older models: MT-XL, MT-X,  
MT6000, MT7, MT7000)  
Reichert Ultracut E, S, and T  
DuPont MT6000

Set-up in minutes

“No Wait” removal from microtome,  
even when cold

“Through-the-Wall” Arm

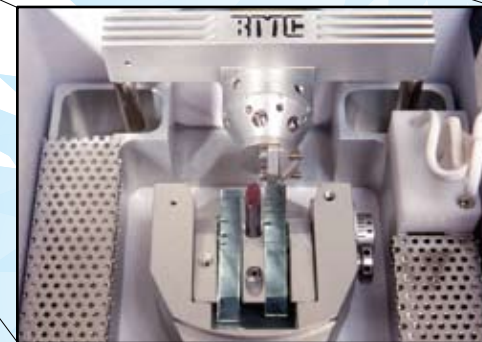

Temperature stability  
± 0.1 C

Temperature range  
+35 to -185 C

Excellent thin  
section quality

Holds two glass or diamond knives  
and trimming tool

Lowest LN<sub>2</sub> consumption,  
9 liters lasts about 11 hours

### Boeckeler®

Boeckeler Instruments, Inc. • 4650 South Butterfield Drive • Tucson, AZ • 85714 • U.S.A.  
Phone (520) 745-0001 • Fax (520) 745-0004 • E-Mail info@boeckeler.com • Web www.rmcproducts.com

© 2008 Boeckeler Instruments, Inc. Boeckeler is a registered trademark of Boeckeler Instruments, Inc., Tucson, Arizona. All specifications subject to change without notice.